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SCHOTT® Instruments
Combination Electrode

CI 60

Chloride Combination Electrode CI 60

**Accuracy when
going to press**

The use of advanced technology and the high quality standard of our products are the result of a continuous development. This may result in differences between this operating manual and your electrode. Also, we cannot guarantee that there are absolutely no errors in this manual. Therefore, we are sure you will understand that we cannot accept any legal claims resulting from the data, figures or descriptions.

**Warranty
declaration**

In case of proper handling, guarantees the quality of the electrode for one year.

Copyright

© Mainz 2012, SI Analytics GmbH,
Reprinting - even as excerpts - is only allowed with the
explicit written authorization of SI Analytics GmbH, Mainz.
Printed in Germany.

View SCHOTT® Instruments Combination Electrode CI 60

Commissioning

1	Remove the protection cap.
2	Pull the closing ring downwards so that the refilling opening of the reference electrode is free.
3	Fill ELY/BR/503 electrolyte solution into the refilling opening.
4	Wipe the shaft using a clean paper towel.
5	Place the combination electrode into deionized water for approx. 10 minutes.
6	Before use, place the combination electrode into diluted standard solution for approx. 5 minutes. In order to achieve a short response time, the concentration of the standard solution should be slightly below the measuring range expected.

The combination electrode is now ready to measure. Always leave the refilling opening of the reference electrode open while measuring.

Calibration

Calibrate according to the data given in the operating manual of the measuring instrument and the analysis specification.

Sample preparation

Add 2 % ISA/FK solution. This sample conditioning solution creates optimum conditions for measuring. It provides a constant ionic strength and similar diffusion potentials at the reference electrode in standard and test sample.

If you would like to have more detailed information concerning sample preparation and measuring procedures, SI Analytics provides a large number of application reports for various applications.

Response times

The response time depends on the concentration range.

It is:

- several seconds with high concentrations,
- several minutes near the detection limit.

The measured value is stable if the value does not change by more than 0.1 mV within 30 seconds.

Interferences

- Complexes with Bi^{3+} , Cd^{2+} , Mn^{2+} , Pb^{2+} , Sn^{2+} , Tl^{3+}
- Reducing agents
- Interfering ions: 10 % error with the following ratios of concentration: (Ratio of concentration = Interfering ion / Measuring ion)

OH^-	Br^-	I^-	S^{2-}	CN^-	NH_3	$\text{S}_2\text{O}_3^{2-}$
80	3×10^{-3}	5×10^{-7}	1×10^{-6}	2×10^{-7}	0.12	0.01

Aging

Please note that every electrode undergoes a natural aging process. The response time increases and the slope decreases with the age of the electrode. The following factors shorten the lifetime considerably:

- Incorrect storing
- Special measuring conditions (e.g. organic solutions, frequent measuring with high concentrations of interfering ions)
- High temperatures
- High changes in temperature

The warranty does not cover damage caused by measuring conditions and mechanical damage.


Maintenance

- For combination electrodes: replenish used reference electrolyte solution.
- In case of increased response times of the electrode perform the following activities:
 - Polish the membrane surface using an SF/K polishing strip. To do so, moisten a polishing strip with deionized water, and with circular movements polish the membrane surface for approx. 30 seconds.
 - Subsequently place the electrode into diluted standard solution for approx. 5 minutes.

Storage

- Between two measurements, place the combination electrode into diluted standard solution.
- Overnight up to one week: Push the closing ring upwards to close the refilling opening. Then rinse the combination electrode with deionized water, blot it dry using a clean paper towel and store it in a dry place.
- For more than a week: Remove the electrolyte solution from the reference electrode and rinse the combination electrode with deionized water, blot it dry using a clean paper towel and put on the protection cap. Store the combination electrode in a dry place.

Typical calibration line of a chloride electrode

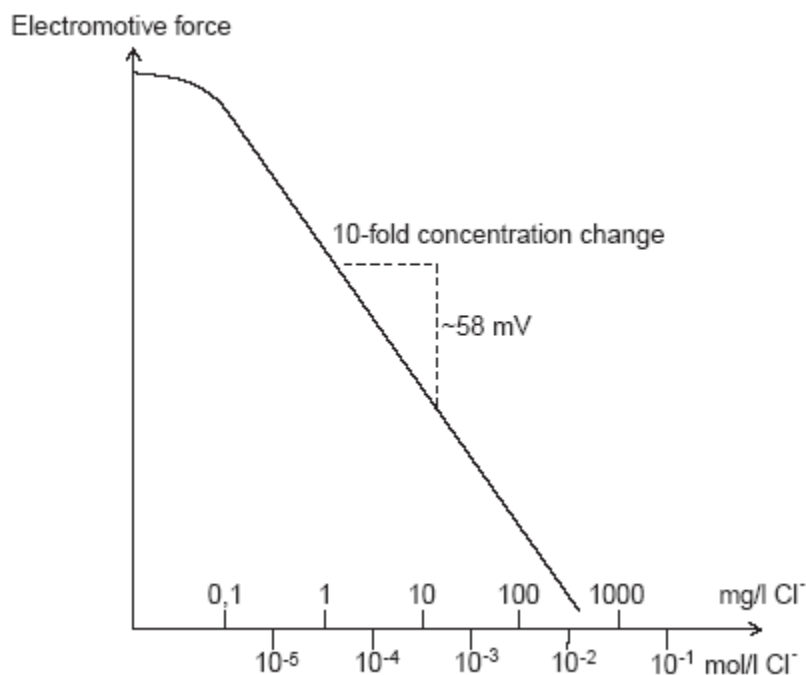Technical data

SCHOTT® Instruments ISE - Combination Electrode Cl 60

Membrane material Pressed piece of ion crystal

Measuring range 2 ... 35000 mg/l Cl⁻ (10⁻⁵ ... 1 mol/l Cl⁻)

Reproducibility ± 2 %

pH range 2 ... 12 (see section INTERFERENCES)

Temperature range 0 ... 80 °C / up to 100 °C for a short time

Membrane resistance < 1 MΩ

Length 153 mm (including 33 mm connection head)

Diameter Shaft: 12 mm
Connection head: 16 mm

Notizen / Notes

SI Analytics

a **xylem** brand

SI Analytics GmbH

Hattenbergstrasse 10

55122 Mainz

Germany

Phone: +49 (0) 6131/66 5111

Telefax: +49 (0) 6131/66 5001

E-mail: support.si-analytics@xyleminc.com

Internet: www.si-analytics.com

T814841